

ClinActiv™ ⊕ MCM™

Next generation of pressure care management with Immersion™ Therapy and Microclimate Management™ for enhanced protection against skin breakdown



Features and benefits

- **Five individual low-pressure zones**, including a dedicated heel zone
- **Real-time pressure adjustment** controlled by the internal Immersion sensor adjusts internal cushion pressures instantly according to patient size, shape, weight or position
- **Flexibility** in choosing between **continuous** or **alternating low-pressure** therapy modes
- **Continuous Microclimate Management** limits the effects of temperature and moisture build-up
- **Dedicated heel zone** with ultra-low pressure for optimal protection of the vulnerable heel area, and with enhanced lower-leg support
- **One touch P-MAX feature** provides fast inflation to for easier patient repositioning and transfer
- **Patient Egress Alarm** informs nursing staff of any undesired patient egress from the bed

ClinActiv ⊕ MCM

Pressure redistributing therapy system


The ClinActiv ⊕ MCM mattress replacement system from Hill-Rom helps to prevent and treat pressure ulcers in low to very high risk adult patients by providing either alternating or continuous low pressure redistribution.

The patented MCM system also helps to manage the patient microclimate and reduces the risk of tissue breakdown as a result of excessive heat and/or moisture.

Choose between Continuous Low Pressure (CLP) or Alternating Low Pressure (ALP) therapy. The Hill-Rom ClinActiv ⊕ MCM provides real immersion therapy in both modes to achieve and maintain the lowest possible interface pressures.

In both modes, the precise level of pressure redistribution is continuously and automatically controlled by the internal immersion sensor which adjusts internal cell pressures according to each patient's size, weight, morphology and position.

Twenty-four-hour pressure ulcer management can be achieved with the addition of the Reflex seat cushion which connects to the same air supply unit, yet can be sealed for greater patient independence whilst out of bed.

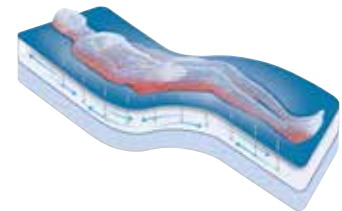

Technical Specifications

Mattress Replacement

Length	203 cm
Width	92 cm
Height	18 cm
Weight	13.9 kg

Control Unit

Length	31 cm
Width	22 cm
Height	28 cm
Weight	4.6 kg

Safe Working Load

250 kg

Recommended therapeutic weight

30 to 150 kg

Cover material

Polyurethane coating on polyester material, 360° zipper, welded seams, low friction, two-way stretch, moisture vapour permeable, bacteriostatic, fungistatic and antimicrobial. Can be wiped and washed.

Hill-Rom is a leading global medical technology company with more than 10,000 employees worldwide. We partner with healthcare providers in more than 100 countries, across all care settings, by focusing on patient care solutions that improve clinical and economic outcomes in five core areas: **Advancing Mobility, Wound Care and Prevention, Patient Monitoring and Diagnostics, Surgical Safety and Efficiency and Respiratory Health.** The people, products and programs of Hill-Rom work towards one mission: **Every day, around the world, we enhance outcomes for patients and their caregivers.**

Not all products/options are available in all countries. For further information about our products and services, please contact your local Hill-Rom representative or visit our webpage:

www.hill-rom.com

The ClinActiv+ MCM surface is intended to help prevent and treat pressure ulcers Stage I, II, III and IV in low to very high risk patients, with a minimum weight of 30 kg and a minimum height of 1.46 m, in a range of healthcare facilities.

Class IIa

Notified Body: LNE/G-Med CE0459

Manufacturer: Hill-Rom S.A.S. – BP 14 ZI du Talhouët 56330 Pluvigner France

This medical device is a regulated health product which, pursuant to such regulation bears a CE mark. Hill-Rom recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices. Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2017 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.

Doc. No: 5EN119301-02, 18 May 2017

Enhancing outcomes for patients and their caregivers:

Hill-Rom