

# Hill-Rom® 900 Elite bed

Improve safety. Simplify tasks.



## Features and benefits

The **lowering of the split siderails** can be activated **using one hand** and feature a smooth, rapid self-lowering mechanism allowing for quick access to patients

Stable support is provided during bed exit via an **ergonomically-placed side egress grip handle**

**Intuitive bed controls** require minimal in-service staff training and help minimise the risk of user's error

**One-Touch Side Egress** positions the bed for safer and easier bed exit

**AutoContour™ and SlideGuard™** minimise patient migration in bed, reducing abdominal compression during bed articulation

**EasyChair™ and flat position function** offers patients the ability to quickly reach a comfortable seated position

**Low height indicator** warns the caregiver if the bed is not in the lowest position, contributing to a safer patient environment

**5th castor** facilitates easier bed transport in any direction

**One-Touch Shock positioning** flattens the upper deck and places the bed into a Trendelenburg position, should patient complications occur

**Intelligent night light** enhances operational safety for both patients and caregivers

Alerts caregivers of potential bed exit via the **3-Mode Bed Exit Alarm**, available with remote monitoring capabilities\*


\* Option available in selected markets

Technical Specifications		Hill-Rom 900 Elite
Safe Working Load	220 kg	●
	250 kg	○
Bed weight		144 kg
Height adjustment		38.9 to 76.8 cm
SlideGuard sliding backrest		up to 12 cm
Protection class		IPX4
Bed angles	Head section	0 to 65°
	Thigh section	0 to 28°
	Foot section	- 22 to -3°
	Trendelenburg/Reverse Trendelenburg	-17 to +17°
Sleep deck	Width	90 cm
	Length	200 cm
Sleep deck surface		High Pressure Laminate (HPL)
Siderail height		39.3 cm
Castors	150 mm Integral	●
	125 mm double band **	○
	150 mm double band	○
Steering castor	At head end	●
	At foot end	○
5th castor		○
Brake and steer pedals	At foot end	●
	At foot and head ends **	○
Brake-off indicator		○
3-Mode Bed Exit Alarm		○
Nurse call		○
Controls	Integrated in siderails	●
	Not integrated in siderails (pendant)	○
	Pendant **	○
	Control on flexible arm **	○
	Caregiver foot end pendant**	○
	Bilateral high-low foot control	○
Head section	Mobile	●
	Stationary	○
Low height indicator **		●
EasyChair and flat position function		●
Intelligent night light		○
One-Touch Shock positioning **		○
Bed extension with linen holder		○

\*\* Only available in combination with selected options.  
For more information please contact your local Hill-Rom representative.

- = not available ● = standard ○ = optional


## Hill-Rom 900 bed with split siderails

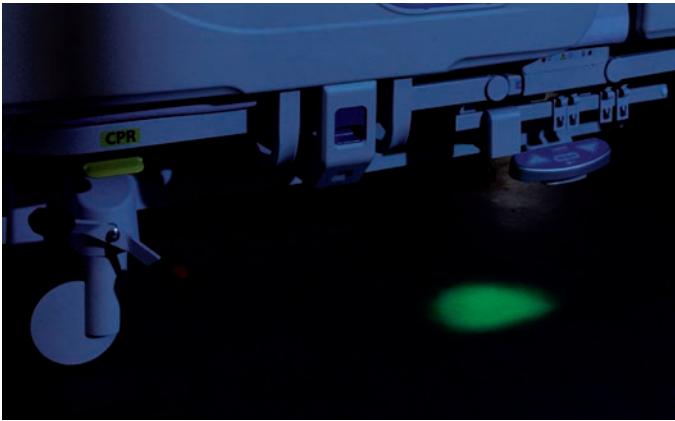

Optional intelligent night light



Low height indicator and One-Touch Side Egress positioning

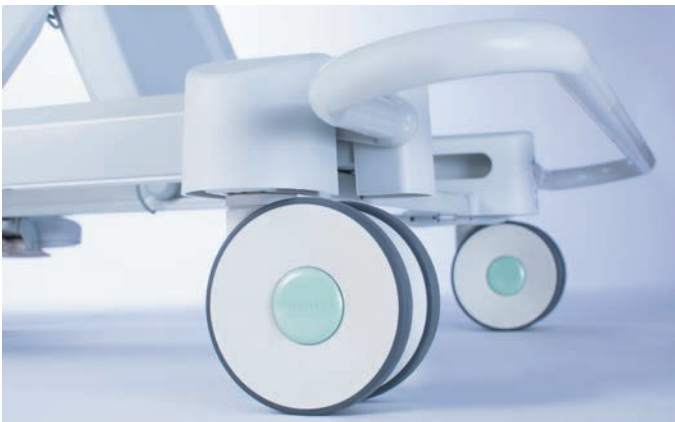

Optional 125 mm double band castors\*\*



3-Mode Bed Exit Alarm



Care grip handles and integrated controls



Optional bilateral high-low foot control with automatic lock-out and 5th steering castor

### Additional standard features

- Bilateral Mechanical CPR lever
- 6 attachment points for patient restraints
- Lifting/IV pole holders
- Bilateral accessory holder
- Electric cord holder
- Adjustable mattress retainers
- Battery
- Removable endboards
- Line-of-site backrest and Trendelenburg indicators
- One-Touch Side Egress positioning
- Central brake
- Service required indicator

Hill-Rom is a leading global medical technology company whose products, services and more than 10,000 employees worldwide help people get better care inside and outside the hospital. Our innovations in five core areas – **Advancing Mobility, Wound Care and Prevention, Patient Monitoring and Diagnostics, Surgical Safety and Efficiency, and Respiratory Health** – improve clinical and economic outcomes and ensure caregivers in more than 100 countries have the products they need to protect patients, speed up recoveries and manage conditions. [Every day, around the world, we enhance outcomes for patients and their caregivers.](#) Learn more at [hill-rom.com](http://hill-rom.com).



Hill-Rom is committed to the concept that mobilising people early and supporting their independence, improves life - in the hospital and at home.

The Hill Rom 900 bed is intended for use in intensive, acute and ambulatory care settings for patients with a weight  $\geq 40$  kg, height  $\geq 146$  cm and BMI  $\geq 17$ .

Class I

Manufacturer: Hill-Rom S.A.S. – BP 14 ZI du Talhouët 56330 Pluvigner France

This medical device is a regulated health product which, pursuant to such regulation bears a CE mark. Hill-Rom recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2018 Hill-Rom Services, Inc. ALL RIGHTS RESERVED. Doc. No: SEN129312-08, 05 December 2018

Not all products/options are available in all countries.  
For further information about our products and services, please contact your local Hill-Rom representative or visit our webpage:

[www.hill-rom.com](http://www.hill-rom.com)

Enhancing outcomes for patients and their caregivers:

**Hill-Rom**